



### AC 1200 - AC 1300 - AC 1400 feeding bottles and animal cages washers



- Independent washing/rinsing system
- Multistage water filtering system, with self-cleaning filter
- AC 1400, single door and pass-through versions available
- High capacity cabinet for detergents
- Easy access to service components

#### Chamber dimensions

**AC 1200**  
W 765 x D 900 x H 815 mm  
W 30.12" x D 35.43" x H 32.09"

**AC 1300**  
W 1550 x D 730 x H 815 mm  
W 61.02" x D 28.74" x H 32.09"


**AC 1400**  
W 1550 x D 900 x H 815 mm  
W 61.02" x D 35.43" x H 32.09"



#### AC 1400 Aquatic Tank Processing

*The specialized version and process to ensure the elimination of any residual biofilm and algae left in aquatic tanks.*

### MITO AC 3500 automated tunnel washer

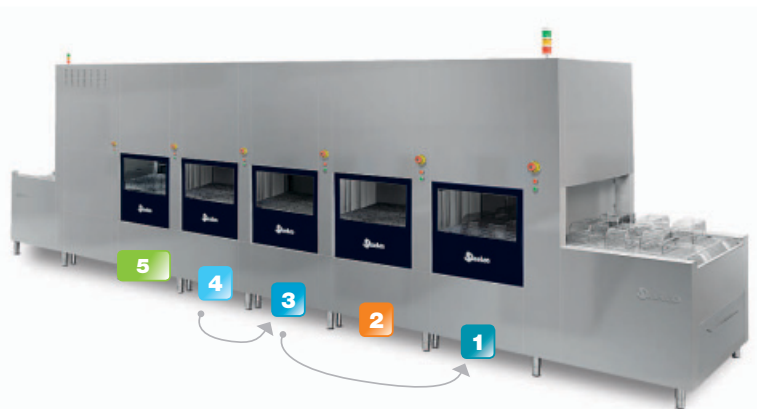

#### Water and energy conservation

Final rinse water (4) is recycled into pre-rinse (3) and from pre-rinse into pre-wash (1)

- Phase cycle:**
- 1) pre-wash
  - 2) chemical wash
  - 3) pre-rinse
  - 4) final fresh hot water rinse
  - 5) hot air drying system

- High productivity
- Multiple chamber modular system for integrated and flexible solutions
- Step or continuous operation versions
- 950mm/35.43" or 1100mm/43.31" belt width version
- Visual inspection of each chamber
- Suitable for integration with robotic technology
- Water filtering system with easy maintenance

